
WelchAllyn

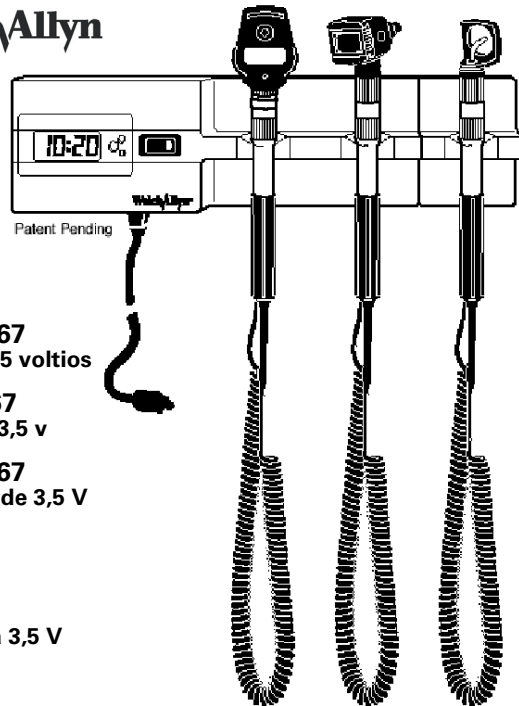Welch Allyn®
767 Diagnostic System
3.5v 767 Wall Transformer

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© 2012 Welch Allyn. All rights are reserved. To support the intended use of the product described in this publication, the purchaser of the product is permitted to copy this publication, for internal distribution only, from the media provided by Welch Allyn.

Welch Allyn assumes no responsibility for any injury to anyone, or for any illegal or improper use of the product, that may result from failure to use this product in accordance with the instructions, cautions, warnings, or statement of intended use published in this manual.

Welch Allyn is a registered trademark of Welch Allyn.

For information about any Welch Allyn product, please call Welch Allyn Technical Support:

USA	+ 1 315 685 4560 800 535 6663	Australia	+ 61 2 9638 3000
Canada	800 561 8797	China	+ 86 216 327 9631
European Call Center	+ 35 3 46 906 7790	France	+ 33 1 55 69 58 49
Germany	+ 49 7477 92 71 86	Japan	+81 42 703 6084
Latin America	+ 1 305 669 9003	Netherlands	+ 31 202 061 360
Singapore	+ 65 6419 8100	South Africa	+ 27 11 777 7555
United Kingdom	+ 44 20 7365 6780	Sweden	+ 46 8 58 53 65 51

Part No. 718764 Rev A 2012-06

Welch Allyn, Inc.
4341 State Street Road
P.O. Box 220 Skaneateles Falls, NY
13153-0220 U.S.A.

www.welchallyn.com

EC	REP
----	-----

Regulatory Affairs Representative
Welch Allyn Limited
Navan Business Park
Dublin Road, Navan
County Meath, Republic of Ireland

TABLE OF CONTENTS	Page
Mounting Instructions	1- 3
Operation	1- 4
Operating Instructions	1- 5
Clock Option.	1- 6
Mounting the Third Handle Module	1- 8
Maintenance.	1- 9
Cleaning and Repair	1- 9
Specifications	1- 10
Guidance and Manufacturer’s Declaration	1- 12

Thank you for purchasing the Welch Allyn 3.5v 767 Wall Transformer. This manual is meant to provide product specifications and instructions for usage and maintenance. The instructions for use should be followed to ensure accurate and reliable service.

The 3.5v 767 Wall Transformer carries a lifetime warranty against all manufacturing defects. **The coiled cords carry a special 10-year warranty against breakage during normal usage.**

Warnings and Cautions

NO MODIFICATION OF THIS EQUIPMENT IS ALLOWED.

USE HOSPITAL-GRADE POWER INLET CORD ONLY.

REPLACE CLOCK BATTERY AS INSTRUCTED WITH Maxell 43 OR EQUIVALENT BATTERY.

GREEN PILOT LIGHT INDICATES THAT THE POWER MAINS IS CONNECTED TO APPLIANCE.

USE THIRD HANDLE MODULE (MODEL #76730) WITH WELCH ALLYN 3.5v 767 WALL TRANSFORMER ONLY.

SERVICE TO BE PERFORMED BY AUTHORIZED WELCH ALLYN REPAIR DEPARTMENT ONLY.

FUSE REPLACEMENT — TYPE T250, 250v SLOW BLOW 0.10 AMPERE (SHOULD BE REPLACED ONLY BY AUTHORIZED WELCH ALLYN REPAIR DEPARTMENT).

TO MINIMIZE THE EXTERNAL HOUSING TEMPERATURE OF THE INSTRUMENT HEAD, ON-TIME MUST NOT EXCEED 2 MINUTES, AND OFF-TIME MUST BE A MINIMUM OF 10 MINUTES.


CAUTION: DISCONNECT THE SYSTEM FROM THE MAINS SUPPLY BEFORE CONNECTING THE THIRD HANDLE (MODEL #76730).


CAUTION: ELECTROMAGNETIC INTERFERENCE RISK. THE DEVICE COMPLIES WITH APPLICABLE DOMESTIC AND INTERNATIONAL STANDARDS FOR ELECTROMAGNETIC INTERFERENCE. THESE STANDARDS ARE INTENDED TO MINIMIZE MEDICAL EQUIPMENT ELECTROMAGNETIC INTERFERENCE. ALTHOUGH THIS DEVICE IS NOT EXPECTED TO PRESENT PROBLEMS TO OTHER COMPLIANT EQUIPMENT OR BE AFFECTED BY OTHER COMPLIANT DEVICES, INTERFERENCE ISSUES STILL MAY OCCUR. AS A PRECAUTION, AVOID USING THE DEVICE IN CLOSE PROXIMITY TO OTHER EQUIPMENT. IN THE EVENT THAT EQUIPMENT INTERFERENCE IS OBSERVED, RELOCATE THE EQUIPMENT AS NECESSARY OR CONSULT MANUFACTURER'S DIRECTIONS FOR USE.

Mounting Instructions

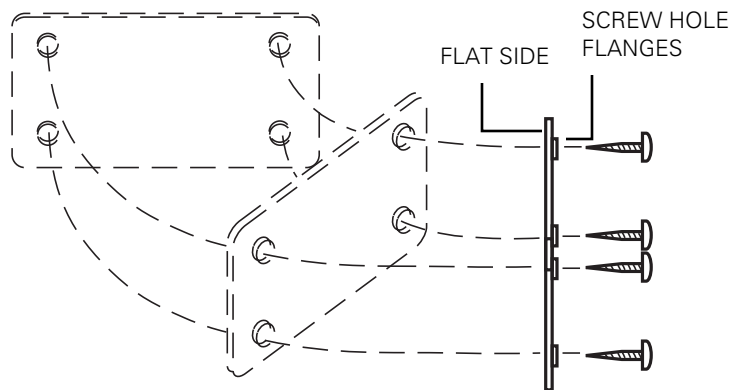
Fasten mounting plate to the wall with enclosed screws, making sure that the screw hole flanges face away from the wall and the smooth surface of the mounting plate is flush with the wall. Depending on the mounting surface, alternative screws may be needed. With the flat side of the plate on the mounting surface, level the plate as shown. Note that the plate should be mounted so that the widest portion of the plate is horizontal (as shown). When mounting screws are tightened, there will be a gap between the screw head and surface of the plate.

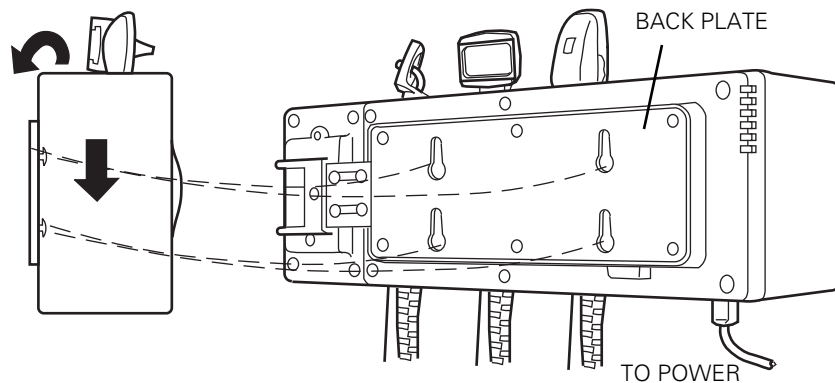Mount the unit by fitting the round holes in the backplate over screw heads on the wall mounting plate, then push downward on the unit to secure it to the wall. This slides the tapered slots in the backplate over the shoulders in the mounting plate. **Plug into 230v, 50 – 60 Hz AC (240v, 50 Hz AC in Australia) and unit is ready for use.**

Operation

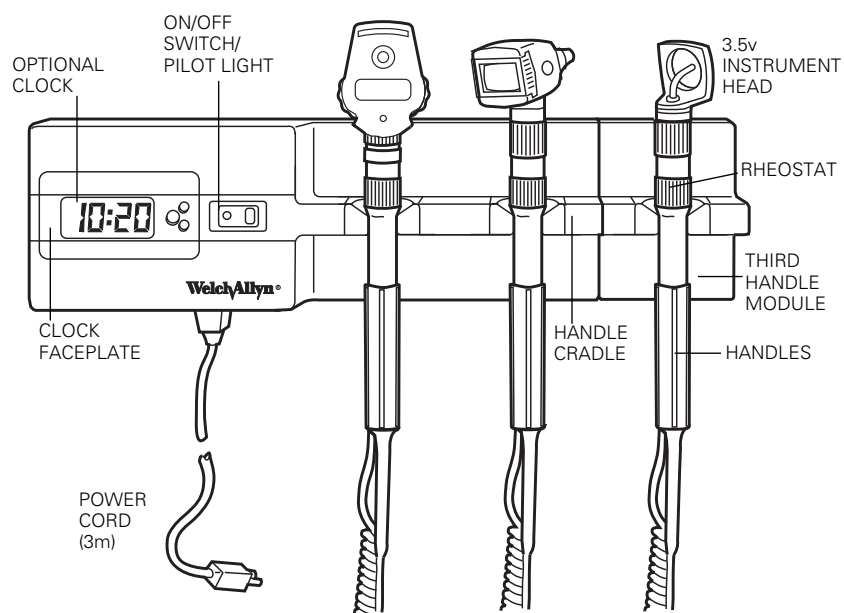Handles Handles will accept any 3.5 v Welch Allyn instrument head.

Rheostat Located on each handle. Turn clockwise to increase light output, turn counterclockwise to decrease light output.

Handle cradle Instrument lights automatically when handle is lifted from the cradle engaging the OptiSense™ Optical Sensor. When handle is returned to the cradle the OptiSense Optical Sensor is disengaged automatically turning the instrument off.

ON-OFF switch/pilot light Indicates the transformer is on (power mains is applied).

Clock Option on this model. Refer to page 1-6 for operating instructions.

Third Handle Module Option on this model. Refer to page 1-8 for mounting and operating instructions.

Power cord Can be plugged into any **230 v, 50 – 60 Hz (240 v, 50 Hz AC in Australia)** outlet.

Operating Instructions


Handles will accept any Welch Allyn 3.5v instrument head.

After plugging power cord into outlet, turn power switch to ON position, illuminating the pilot light.


With power switch in ON position, instruments will automatically light when handles are lifted from the handle cradle and the OptiSense Optical Sensor is engaged. Adjust rheostat on handle until desired light output is obtained. Turning rheostat clockwise increases light output, counterclockwise decreases light output. Note: Rheostat does not turn instrument completely off. Returning handle to cradle will disengage the OptiSense Optical Sensor and turn the instrument completely off.

A built-in, automatic voltage regulator provides maximum illumination and lamp life.

Instruments will automatically turn off when the handle is placed back in the cradle. The Welch Allyn 767 Wall Transformer can be turned off by simply turning the power switch to the OFF position. Pilot light alone draws current when instruments are not in use. When unit is not in use for more than a few hours, turn power switch to OFF position to ensure longer operating life. Unit can remain plugged in permanently. The clock is powered by a battery and will keep time when power switch is turned off.

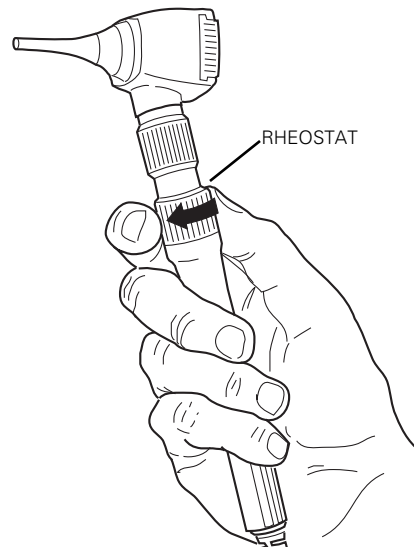Clock Option

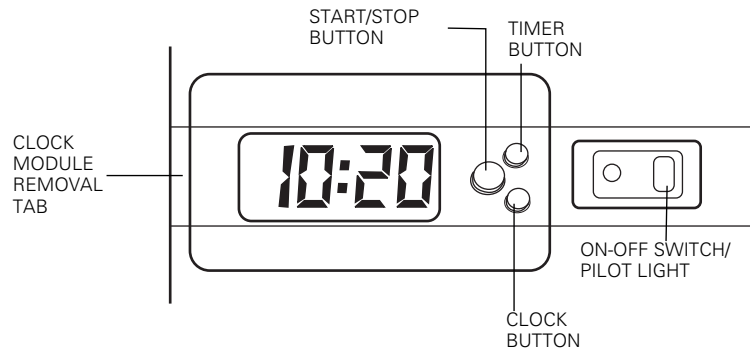Start/Stop button Use this button to START AND STOP seconds count. Press once to start the seconds. Press again to stop seconds.

Clock button Use this button to return from seconds mode to hour/minute by pressing once.


Timer button Use this button to change to seconds mode (timer).

Clock module removal tab Pull forward to remove clock module and gain access for changing battery.

The clock, powered by a battery, operates separately from the main power supply. For this reason the time will be displayed when the main unit is turned OFF or unplugged from the wall power socket.

Engage Seconds Mode:


Press Timer button to engage. Press Start/Stop button to begin seconds mode. Press the Clock button to return the Hour/Minute display.

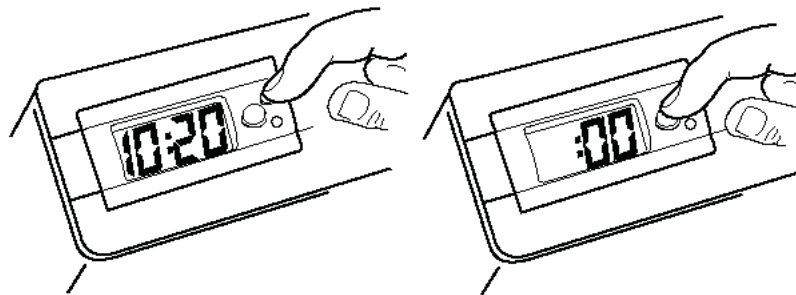
Change Time:

The following steps should be followed when resetting the time:

1. Remove clock module by prying the clock faceplate from unit housing.

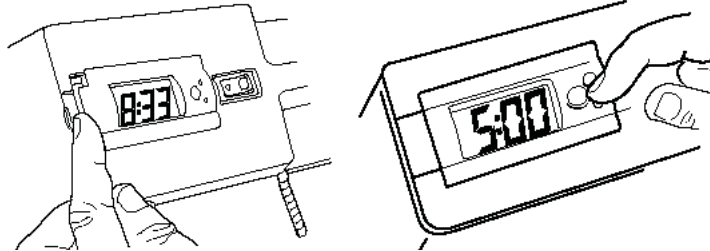
2. With clock in clock mode, press Start/Stop button for 5 seconds until time on clock begins to flash.
3. Using Hour Set and Minute Set buttons on back of clock module, set time to correct values.
4. Press Start/Stop button again to return clock to regular operation.

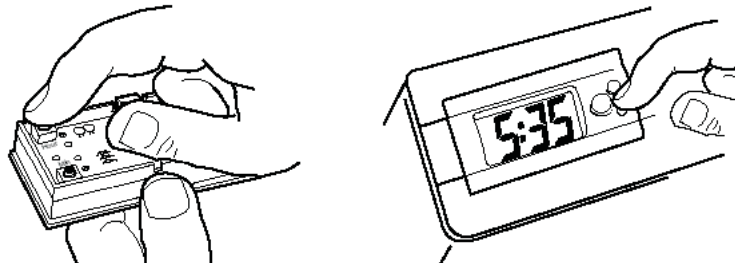5. Replace clock module into unit by lining up edges and pressing into main housing.

Change Clock Battery:

1. Remove clock module by prying the clock faceplate from unit housing.
2. Remove battery by pushing battery out of retaining sleeve. For best results, use a pen or pencil to pry against the faceplate housing wall.
3. Insert new battery. Welch Allyn recommends Maxell 43 or equivalent.

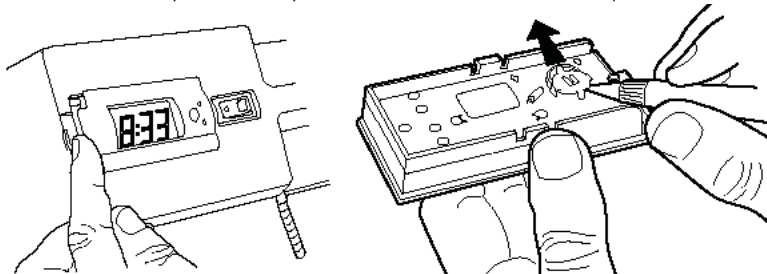
Mounting the Third Handle Module

A Phillips screwdriver is needed to connect the Third Handle Module.


Mounting the Third Handle Module:

1. Disconnect main unit power supply.
2. Remove main unit from the wall by lifting up and out (instructions found on page 1-3 of this manual).
3. Attach stabilization bar to back of Third Handle Module (as shown using two enclosed star washers and screws).

Step 3:

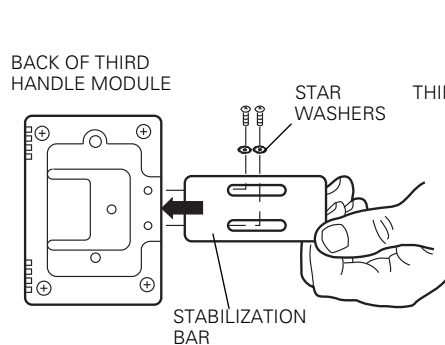
Step 5:

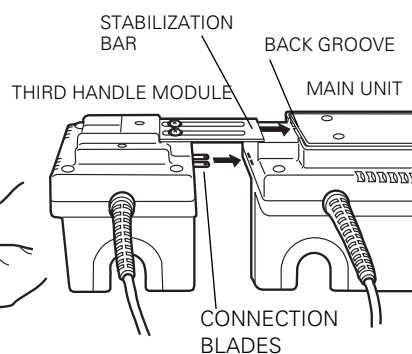4. Line up the metal stabilization bar with the groove on the back of the main unit. Also, line up the connecting blades with the blade openings on the side of the main unit housing.
5. Slide Third Handle Module into the main unit until firmly in place.

STEP 6:

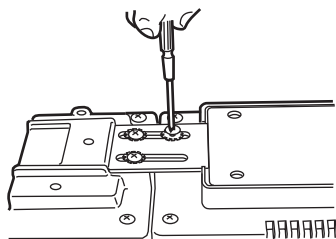
6. Secure Third Handle Module in place (as shown using two enclosed star washers and screws).
7. Place unit on the wall following the instructions found on page 1-3 of this manual and reconnect main unit power supply.
8. Place instrument on third handle and lift from cradle. This will test for correct attachment.

Maintenance

Like any other piece of electrical equipment, periodic electrical inspections should be performed by qualified personnel. Welch Allyn recommends inspections every six months, more often if used under adverse conditions.

Cleaning and Repair

Do not sterilize. May be cleaned by wiping with a dry cloth. Avoid using abrasive cleaning materials on clock faceplate.

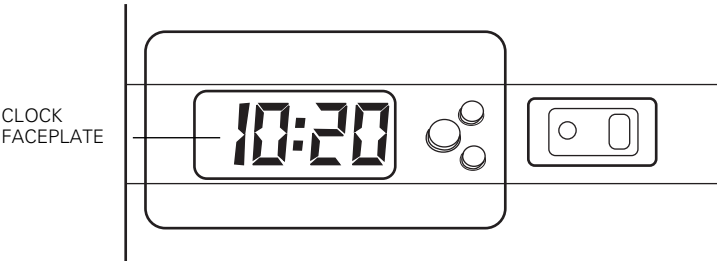REPAIR:

USA + 1 315 685 4560 800 535 6663	Australia + 61 29 638 3000
Canada 800 561 8797	China + 86 216 327 9631
European Call Center + 353 46 906 7790	France + 33 15 569 5849
Germany + 49 747 792 7186	Japan +81 42 703 6084
Latin America + 1 305 669 9003	Netherlands + 31 202 061 360
Singapore + 65 6419 8100	South Africa + 27 11 777 7555
United Kingdom + 44 207 365 6780	Sweden + 46 85 853 6551

Specifications

Model No.

76712	3.5 v 767 Wall Transformer EUR (two-handle model without clock, instrument heads not included),
76714	3.5 v 767 Wall Transformer UK (two-handle model without clock, instrument heads not included),
76716	3.5 v 767 Wall Transformer AUS (two-handle model without clock, instrument heads not included),
76722	3.5 v 767 Wall Transformer EUR (two-handle model with clock, instrument heads not included),
76724	3.5 v 767 Wall Transformer UK (two-handle model with clock, instrument heads not included),
76726	3.5 v 767 Wall Transformer AUS (two-handle model with clock, instrument heads not included),
76730	3.5 v 767 Wall Transformer Third Handle Module (instrument heads not included).

Applied Parts: The handles (with heads attached) are considered to be applied parts.

Power Cord:

Input:	230 v, 50 – 60Hz. (0.5A) EUR/UK
	240 v, 50 – 60Hz. (0.5A) AUS

Supplied with two coiled cords, two handles, and power cord.

Classification: Class I

Leakage current is less than 10 microamps from any exposed metal part.

Fuse: Replace with fuse specified on page 1-1

Weight: Weight with mounting plate is 4.5 pounds (2.04kg). (NOTE: This is weight without instrument heads.)

Dimensions: H=4" (10.16cm) D=4" (10.16cm) W=12" (30.48cm)

Approvals IEC 60601-1, UL 60601-1, IEC 60601-1-2, CAN/CSA C22.2 No 601.1-M90, AS/NZS3200

Type BF Applied Parts

Caution In this manual, indicates conditions that could damage equipment or other property.

Warning In this manual, indicates conditions that could lead to illness, injury, or death.


The CE mark on this product indicates it has been tested to and conforms with the provisions noted within the 93/42/EEC Medical Device Directive.


Regulatory Affairs Representative
Welch Allyn Limited
Navan Business Park
Dublin Road, Navan
County Meath, Republic of Ireland

Operating:

Temperature Range Humidity Range

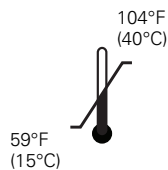
Atmospheric Pressure Range

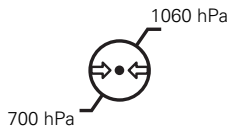
Transport/storage:

Temperature Range Humidity Range

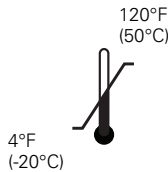Atmospheric Pressure Range

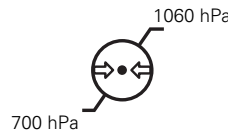Guidance and Manufacturer's Declaration

EMC Compliance

Special precautions concerning electromagnetic compatibility (EMC) must be taken for all medical electrical equipment. This device complies with IEC 60601-1-2.

- All medical equipment must be installed and put into service in accordance with the EMC information provided in this document.
- Portable and mobile RF communications equipment can affect the behavior of medical electrical equipment.

This device complies with all applicable and required standards for electromagnetic interference.

- It does not normally affect nearby electronic equipment and devices.
- It is not normally affected by nearby electronic equipment and devices.
- It is not safe to operate this device in the presence of high-frequency surgical equipment.

It is a good practice to avoid using the device in extremely close proximity to other electronic equipment.

Emissions and Immunity Information

Electromagnetic Emissions

The device is intended for use in the electromagnetic environment specified below. The customer or user of the device should assure that it is used in such an environment.

Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	Group 1	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

Electromagnetic Immunity

The device is intended for use in the electromagnetic environment specified below. The customer or user of the device should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV line(s) to line(s) ±2 kV line(s) to earth	±1 kV line(s) to line(s) ±2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions, and voltage variations on power supply input lines. IEC 61000-4-11	>95% dip for 0.5 cycle 60% dip for 5 cycles 30% dip for 25 cycles >95% dip for 5 seconds	>95% dip for 0.5 cycle 60% dip for 5 cycles 30% dip for 25 cycles >95% dip for 5 seconds	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or battery.
Power frequency (50 – 60Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

Electromagnetic Immunity

The device is intended for use in the electromagnetic environment specified below. The customer or user of the device should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms 150 kHz to 80 MHz	<p>Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance</p> $d = (1.2) \sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m 80 MHz to 2.5 GHz	<p>$d = (1.2) \sqrt{P}$ 80 MHz to 800 MHz</p> <p>$d = (2.3) \sqrt{P}$ 800 MHz to 2.5 GHz</p> <p>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 

Note 1: At 80 MHz and 800 MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the device.

^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Recommended Separation Distances Between Portable and Mobile RF Communications Equipment and the device

The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or user of the device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the device as recommended below, according to the maximum output power of the communications equipment.

Rated Max. Output Power of Transmitter (W)	Separation Distance According to Frequency of Transmitter (m)		
	150 kHz to 80 MHz $d = (1.2) \sqrt{P}$	80 MHz to 800 MHz $d = (1.2) \sqrt{P}$	800 MHz to 2.5 GHz $d = (2.3) \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

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Printed in U.S.A.